

**flexicare**<sup>®</sup>  
TOTAL QUALITY - TOTAL CARE

feel  
confident

Continence Care

With over 30 years of experience in design, development and manufacturing, Flexicare has developed a world class reputation.



“Our commitment is to deliver on our promises”



**Customer Care** Our Customer Care Department is available to respond to your enquiries and to discuss your order requirements. Our experienced Sales team can offer advice and educational support within the rapidly changing hospital environment. We also provide online information regarding our product range.

**Training** We are dedicated to supporting Hospitals in terms of educational support and training of all medical staff with understanding of the principles underlying our product range. We work closely with a number of University Teaching Hospitals and Clinicians to develop full training packages.

**Design and Manufacture** With our experienced in-house Research and Design team we are able to ensure that new products are developed to meet our customers' changing needs. All Flexicare products are rigorously tested to ensure they are of the quality and standard that our customers expect.

**Accreditation** Flexicare is an ISO 13485 Registered Manufacturer of sterile medical devices and surgical products with CE mark accreditation.

FDA Registration No: 8021995



**International Distributors** Flexicare distributes to over 72 countries worldwide through a network of distributors and its own subsidiaries. Flexicare has offices across Europe, the Middle East, Far East and the USA which are able to respond to customers' needs, providing a fast and effective service.

Dedicated to **supporting hospitals** providing comprehensive training

Flexicare distributes to **over 72 countries** worldwide

## CONTENTS

Flexicare offers a wide range of continence care products to meet different patients' needs so that they can continue with their normal daily activities. All Flexicare products are designed and manufactured with every user's requirements taken into consideration, achieving a high level of user confidence and comfort.

### Catheters

FlexiCath All Silicone Foley Catheter \_\_\_\_\_ 4

### Drainage Bags

Flexicare Leg Bags \_\_\_\_\_ 6

Discreet Leg Bags \_\_\_\_\_ 7

Non-Drainable and Drainable 2L Bags \_\_\_\_\_ 8

F2-EZ Once Drainable 2L Bag \_\_\_\_\_ 9

F2-EZB Once Drainable 2L Bag \_\_\_\_\_ 10

Urimeter \_\_\_\_\_ 11

Post Operative Bag \_\_\_\_\_ 12

### Accessories

FlexiHang Bag Hanger \_\_\_\_\_ 13

EZ-Flow Catheter Valve \_\_\_\_\_ 13

Bag Stand \_\_\_\_\_ 14

FlexiSleeve Leg Bag Holder/Sleeve \_\_\_\_\_ 14

Leg Straps \_\_\_\_\_ 14

Home Packs \_\_\_\_\_ 15

Contact Us \_\_\_\_\_ 15



FlexiCath is a 100% silicone 2-way biocompatible indwelling foley catheter designed to eliminate the risk of allergic reactions and provide maximum long-term patient comfort.

# flexiCath

### Smooth Rounded Tip

Allows for easy insertion and provides maximum patient compliance.

### Large Twin Catheter Eyes

The smooth moulded drainage eyes reduce the risk of snagging, ensuring easier insertion when in place. The large eyes allow quick fluid drainage and reduce risk of blockages.

### Large Inner Lumen

FlexiCath has a large inner lumen with a round, smooth surface reducing the risk of bacterial colonization and encrustation, allowing quick drainage of fluids and any possible clots.

### Smooth Valve Surface

The valve sleeve provides a smooth surface against the patient's skin reducing the risk of irritation.

### Radiopaque Line

The integral blue line running the length of the catheter is radiopaque, allowing clinicians to see the position of the catheter under x-ray.



### Spherical Retention Balloon

The ribbed design of the balloon gives perfect symmetry providing a good position and seal when in situ. The silicone balloon has the ability to contract back to its unfilled streamline shape providing easier removal even after 30 days of use.

### Flexible Shaft

Designed to provide optimum stiffness for easy insertion whilst maintaining true flexibility, giving users freedom of movement without discomfort.

### 100% Silicone


FlexiCath is made of 100% silicone which minimizes bacterial colonization and encrustation whilst in situ, facilitating easier removal.

### Ribbed Outlet Connector

Provides a secure attachment to a catheter valve or leg bag, minimizing the risk of slippage, disconnection or leakage. Clearly marked with sizing information.

### Recessed Inflation Valve

Prevents accidental emptying of the catheter balloon.

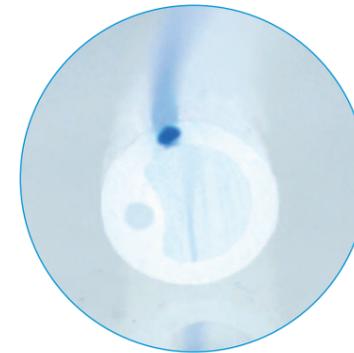

Large inner lumen



Atraumatic tip & large eyes



Ribbed outlet connector

### Paediatric Catheters



Product Code	Color	Ch/Fr	Volume	Quantity
00-410603	Brown	6 Ch/Fr	3-5ml	10
00-410803	Black	8 Ch/Fr	3-5ml	10
00-411003	Grey	10 Ch/Fr	3-5ml	10

### Female Catheters



Product Code	Color	Ch/Fr	Volume	Quantity
00-421205	White	12 Ch/Fr	10ml	10
00-421405	Green	14 Ch/Fr	10ml	10
00-421605	Orange	16 Ch/Fr	10ml	10

### Male Catheters



Product Code	Color	Ch/Fr	Volume	Quantity
00-411205	White	12 Ch/Fr	10ml	10
00-411405	Green	14 Ch/Fr	10ml	10
00-411605	Orange	16 Ch/Fr	10ml	10
00-411630	Orange	16 Ch/Fr	30ml	10
00-411805	Red	18 Ch/Fr	10ml	10
00-411830	Red	18 Ch/Fr	30ml	10
00-412005	Yellow	20 Ch/Fr	10ml	10
00-412030	Yellow	20 Ch/Fr	30ml	10
00-412205	Purple	22 Ch/Fr	10ml	10
00-412230	Purple	22 Ch/Fr	30ml	10
00-412405	Blue	24 Ch/Fr	10ml	10
00-412430	Blue	24 Ch/Fr	30ml	10
00-412605	Pink	26 Ch/Fr	10ml	10
00-412630	Pink	26 Ch/Fr	30ml	10

## Leg Bags - Flexicare Range

Flexicare leg bags are available in 350ml, 500ml and 750ml capacities, meeting patient requirements. All three leg bag capacities are available in short or long tube lengths and for those requiring extra length, Flexicare offer an extra-long adjustable length tube.

### Easy to Operate Leak-Free Taps

The lever taps fitted to Flexicare bags are 100% leak tested and are easy to operate.

### Needle-Free Sample Port

The needle-free sample port eliminates the risk of needle stick injury, whilst ensuring the bag is leak-free.

### Skin-Friendly Fabric Backing

All of our community leg bags have a durable, non-woven fabric backing to reduce skin irritation and perspiration.

### Secure and Comfortable Fit

Anti-slip, latex-free Lycra Velcro straps provide secure, comfortable support and help maintain the inlet tube and non-return valve in position, preventing kinking of the tube or folding of the NRV.

### Universal Secure Overnight Connection

All outlet taps on Flexicare leg bags are fitted with a silicone overnight bag connection to attach securely and reliably to a bedside night drainage bag.

### Range of Tube Lengths

The table below shows the different tube lengths available.

Tube Description	Short Tube	Long Tube	Extra Long Adjustable Tube
Tube Length	7cm / 2 ¾ inches	30cm / 11¾ inches	60cm / 23 ¾ inches

### Flexicare Leg Bags - Community packs with fabric backing

00-1352	350ml Leg Bag Short Tube - sterile	10
00-2352	350ml Leg Bag Long Tube - sterile	10
00-1502	500ml Leg Bag Short Tube - sterile	10
00-2502	500ml Leg Bag Long Tube - sterile	10
00-3502	500ml Leg Bag Extra Long Adjustable Tube - sterile	10
00-3501	500ml Leg Bag Extra Long Adjustable Tube - non-sterile	10
00-1752	750ml Leg Bag Short Tube - sterile	10
00-2752	750ml Leg Bag Long Tube - sterile	10

### Flexicare Leg Bags - Hospital packs without fabric backing

00-1503L	500ml Leg Bag Short Tube - sterile	100
00-2503L	500ml Leg Bag Long Tube - sterile	100



## Leg Bags - Discreet Range

Flexicare offer a discreet range of leg bags which are contoured to fit the anatomy better, have a neutral colour and text that is harder to see through clothing. This allows the user to wear a wider range of clothes without compromising on lifestyle choices.

# Discreet

### Neutral Colour

Using neutral colours for the backing and print on the leg bags allows users to wear the bag with more confidence by reducing the possibility of printed text or bright colour backing showing through clothing.

### Contoured Fit

The Discreet range of leg bags are designed to fit better and mould to either the calf or thigh, distributing the weight and volume of the urine evenly to minimise visibility of a leg bag in-situ.

### Easy to Operate Leak Free Taps

The lever taps fitted to Flexicare bags are 100% leak tested and are easy to operate.

### Needle-Free Sample Port

The needle-free sample port eliminates the risk of needle stick injury, whilst ensuring the bag is leak free.



The Discreet range of leg bags are available in 350ml, 500ml and 750ml capacity with short or long tube and for those requiring extra flexibility, Flexicare offers an adjustable length tube in the 500ml capacity.

### Range of Tube Lengths

The table below shows the different tube lengths available.

Tube Description	Short Tube	Long Tube	Extra Long Adjustable Tube
Tube Length	7cm / 2 ¾ inches	30cm / 11¾ inches	60cm / 23 ¾ inches

### Discreet Leg Bags - Community packs with fabric backing

00-5352	350ml Discreet Leg Bag Short Tube - sterile	10
00-6352	350ml Discreet Leg Bag Long Tube - sterile	10
00-5502	500ml Discreet Leg Bag Short Tube - sterile	10
00-6502	500ml Discreet Leg Bag Long Tube - sterile	10
00-6504	500ml Discreet Leg Bag Extra Long Adjustable Tube - sterile	10
00-6505	500ml Discreet Leg Bag Extra Long Adjustable Tube - non-sterile	10
00-5752	750ml Discreet Leg Bag Short Tube - sterile	10
00-6752	750ml Leg Bag Long Tube - sterile	10

### Discreet Leg Bags - Hospital packs without fabric backing

00-5503	500ml Discreet Leg Bag Short Tube - sterile	100
00-6503	500ml Discreet Leg Bag Long Tube - sterile	100

## Non-Drainable and Drainable 2L Bags

The requirements of all users are taken into consideration in the design and manufacture of Flexicare urine bags. Achieving a high level of user confidence and comfort is a fundamental criterion in the making of Flexicare's range of urology products. All Flexicare bags are made and tested to meet international standards (BS EN ISO 8669).

# f2 & f4

### Safe and Secure Drainage

The unique design and manufacture of our non-return valves (NRV) allows low pressure drainage without compromising its anti-reflux property, a balance required to allow free flow of urine into the bag and at the same time protect the patient against potential ascending infection.

The stepped connector glides easily into both catheter and sheath providing a secure connection that will withstand the most vigorous user activity.

### Needle-Free Sample Port

The needle-free sample port eliminates the risk of needle stick injury, whilst ensuring the bag is leak free.

### Free and Rapid Flow


The wide bore drainage tube allows unimpaired flow of urine into the bag, optimising drainage and providing observation of urine flow. The tube is both firm and flexible affording comfort to the user.

### Easy to Operate Leak-Free Taps

Both the lever tap and T-tap fitted to Flexicare F4 bags are 100% leak tested and are easy to operate. The lever tap requires a lower force to open and close the tap and this design is often preferred by the elderly or patients with impaired dexterity.



F4L

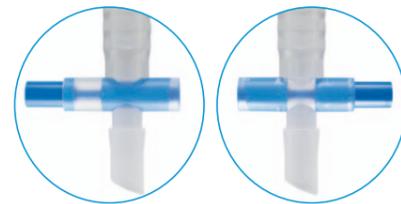

F4



### Overnight Drainage Bags - Community Packs

00-1200	F4 Community Pack - 2L drainable sterile bag with inlet connector, needle-free sampling port and T-tap	10
00-1201	F4L Community Pack - 2L drainable sterile bag with inlet connector, needle-free sampling port and lever tap	10
00-2202c	F2 Community Pack - non-sterile bag, non-drainable with non-return valve and tear to open	10

### Overnight Drainage Bags - Hospital Packs

00-1204	F4 Hospital Pack - 2L drainable sterile bag with inlet connector, needle-free sampling port and T-tap	100
00-1204L	F4 Hospital Pack - 2L drainable sterile bag with inlet connector, needle-free sampling port and lever tap	100
00-2201	F2 Non-Sterile Hospital Pack - non-sterile bag, non-drainable, non-return valve and tear to open	250
00-2202	F2 Sterile Hospital Pack - sterile bag, non-drainable, non-return valve and tear to open	100
00-2206	F3 Non Sterile Hospital Pack - non-sterile bag with non-return valve and T-tap	100

## Once Drainable 2L Bags

F2-EZ is an innovative 2L drainage bag designed to be opened quickly and easily while maintaining controlled drainage of fluids. The twist-snap tap is located at the top of the bag thus eliminating the possibility of splash, spray and cross contamination associated with tear bags. The bag is designed for single use and, once the twist-snap tap has been broken and the bag emptied, it can be easily disposed of.

# f2-EZ



### F2-EZ Twist-Snap Tap

Innovative easy-to-open angled connector with wide bore for rapid drainage.

### Non-Return Valve

Eliminates reflux and maintains directional flow of fluids.

### Opposing Diagonal Eyelet

The positioning of the eyelet allows the bag to be hung for emptying.

### Stepped Connector

This provides a secure connection to all bags and catheters.

00-2203C	F2-EZ Community Pack - non-sterile bag, once-drainable with non-return valve and twist-snap tap	10
00-2203	F2-EZ Hospital Pack - non-sterile bag, once-drainable with non-return valve and twist-snap tap	100

## Once Drainable 2L Bags

F2-EZB features the same specification as F2-EZ but with the twist-snap tap at the bottom of the bag. The bag is designed for single use and, once the twist-snap tap has been broken and the bag emptied, it can be easily disposed of.

# f2-EZB

**Non-Return Valve**  
Eliminates reflux and maintains directional flow of fluids.

**Stepped Connector**  
Provides a secure connection to all bags and catheters.

**Twist-Snap Tap**  
Held in the upright position the twist-snap tap is twisted to break the seal.



**Bag Stand**  
F2-EZB is pictured with the Flexicare Bag Stand. For more information see page 14.

**F2-EZB Twist-Snap Tap**  
Innovative easy-to-open connector with wide bore for rapid drainage.



**Easily Drainable**  
After discarding the tip, the fluid can be easily drained via the outlet.

00-2208C	F2-EZB Community Pack - non-sterile bag, once-drainable with non-return valve and bottom-located twist-snap tap	10
00-2208	F2-EZB Hospital Pack- non-sterile bag, once-drainable with non-return valve and bottom-located twist-snap tap	100

## Urimeter

Urine meters are primarily used with patients post-operatively or those critically ill, where precise monitoring of urine output is necessary. Flexicare have developed a sterile closed collection Urimeter with an integral measuring chamber directly above the in-line collection bag.


**Secure Fastening**  
The integral hangers and straps ensure the Flexicare Urimeter is securely attached to any bed frame.

**Vertical Drainage**  
Vertical drainage ensures complete emptying of the measuring chamber

**Safe and Secure Drainage**  
The unique design and manufacture of our non-return valves (NRV) allows low pressure drainage without compromising its anti-reflux property.

**Hydrophobic Filters**  
The hydrophobic filter allows for the air and any gas build up in the bag to escape without allowing proteins or fluids to pass through.

**Protective Tap Holder**  
The lever tap can be housed in the blue tap holder when not needed to prevent the tap touching the floor or being accidentally opened.

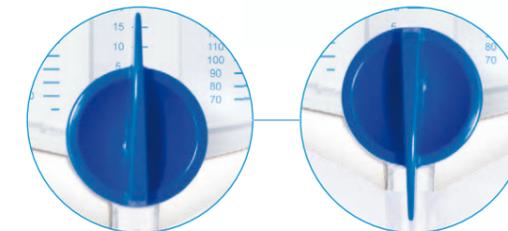

**Easy Emptying**  
The 180° release tap enables the clinician to empty the measuring chamber quickly and without the need to invert the chamber, reducing the risk of contamination.

Urine flows uninterrupted into the main measuring chamber. This chamber is divided into 3 sections with the central chamber measured in clearly defined 5ml increments allowing for accurate measurements up to 65ml. Once full, urine then flows into the second and third chambers which are measured in 10ml increments with a total capacity of 500ml. Overflow urine drains into the collection bag.

The urine can be evacuated from the measuring chamber into the collection bag via an easy to open rotating tap. This in-line system minimizes the risk of contamination and ensures complete emptying of the measuring chamber should further monitoring be necessary.

00-1225	Flexicare Urimeter. 500ml Chamber and 2000ml Bag	20
---------	--	----



## Post Operative Bag

Primarily used with patients who have undergone an operation and are fitted with a urinary catheter. The Flexicare Post Operative Bag is designed for patients with a weakened immune system who require a sterile drainage bag with maximum protection against the risk of ascending infection.

### Integral Bag Hanger

The integral bag hanger can fit securely to any bed or stand. The inlet tube is also maintained in its primary position to allow for direct flow into the drip chamber.



### Drip Chamber

The drip chamber is fitted with a non-return valve and prevents contact between the urine inside the bag and the tube, which eliminates the risk of ascending infection.

### Hydrophobic Filter

The hydrophobic filter allows for any air or gas build up in the bag to escape without allowing proteins or fluids to pass through. The function of the hydrophobic filter also assists in proper filling of the bag without pressure build up.

### Non-Return Valve

The bag is fitted with a second NRV preventing the contents of the bag returning to the drip chamber, eliminating reflux and provides directional flow of fluids.



### Secure Lever Tap

The lever tap can be placed in a protective tap holder which keeps the outlet from possible cross contamination from the floor whilst also protecting against accidental opening.

The bag is connected to the catheter via a needle free sample port connector which provides a secure fit to any foley catheter. The wide bore tubing allows for quick and easy urine drainage, while remaining flexible and kink resistant.

00-1208


Flexicare Post Operative Bag - Sterile. 2000ml Capacity

20

## Accessories

### flexiHang

The FlexiHang bag hanger is a uniquely designed universal 2 litre urine bag hanger for use on beds fitted with or without side bars. FlexiHang can also be secured between the mattress and bed frame in a divan bed. The ridges incorporated into the hook section of the hanger will provide a positive grip keeping the hanger and bag in position. FlexiHang holds the urine bag firmly in place with two snap-fit locks and secures the outlet tube in a vertical position to ensure good flow to the bag.

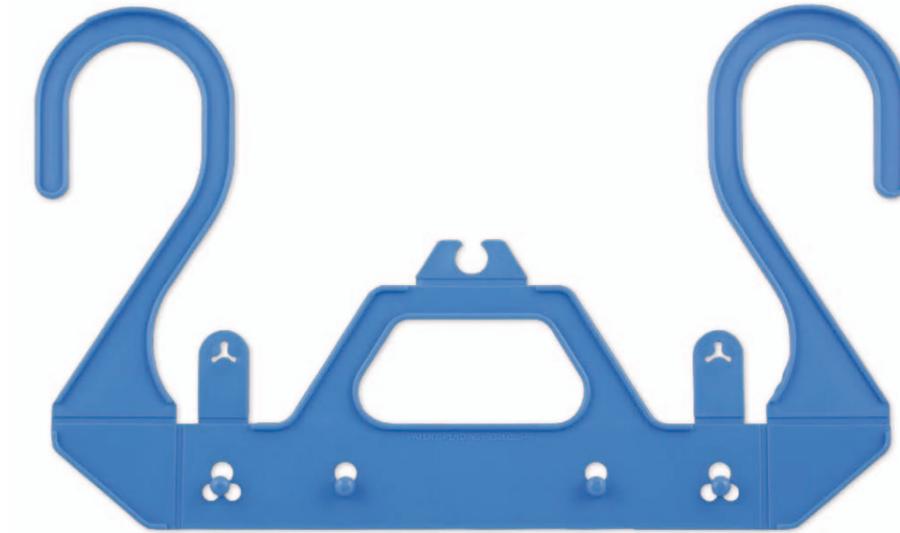

00-0012 FlexiHang Bag Hanger

10

00-0014 FlexiHang Bag Hanger

100

### EZ-flow

The Flexicare EZ-Flow is a discreet lightweight catheter valve designed to offer patients flexibility of choice in urine drainage management.

**100% tested** - All EZ-Flow catheter valves are 100% tested to guarantee leak-free valves.

**Patient Comfort** - The outlet sleeve of the catheter valve is made of soft PVC that will not cause pressure marks or discomfort for the patient.

**Secure Fitting** - the ridged connector ensures a secure fit with the catheter.



00-0060

EZ Flow Catheter Valve

5


## Accessories

### Bag Stand

The Flexicare Bag Stand is tough, durable and designed to hold any 2L bag. The metal stand is plastic coated making it easy to clean and ensuring that there are no sharp edges.



Easy Clean Coating

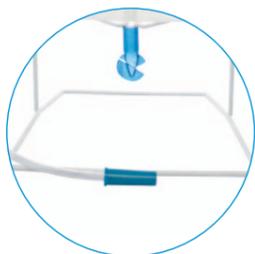

Stable Base



00-0022

Bedside Bag Stand

50

### FlexiSleeve

FlexiSleeve is a knitted leg bag holder to replace leg straps. FlexiSleeve offers more comfort to the bag wearer by distributing the weight of the bag evenly around the leg or thigh. This helps to reduce the risk of pressure sores occurring and FlexiSleeve also helps to dampen the noise of the fluid in the bag.



Increased Comfort



Knitted Weave



00-0035

Small size leg bag holder

3

00-0036

Medium size leg bag holder

3

00-0037

Large size leg bag holder

3


00-0038

Extra-large size leg bag holder

3

### Leg Straps

Flexicare leg straps are made of Lycra which keeps its elasticity over the strap's lifetime. The leg straps feature an anti-slip silicone wave which grips the leg and ensures a secure and comfortable fit.

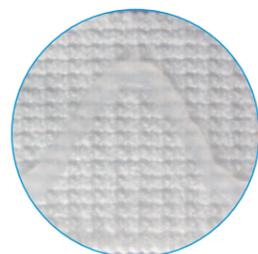

Anti-Slip Silicone Wave



Secure Fastening



00-0032C

Lycra anti-slip leg straps

10

00-0032

Lycra anti-slip leg straps

100

## Home Packs

Flexicare offers home packs to patients leaving hospital and returning home. Several different home packs are available to suit different patient requirements, ensuring that patients have all the items they need during their first week away from hospital.

All home packs include a user booklet, written by a Continence Advisor, designed specifically for patients to ensure they use the products correctly as well as emphasising the importance of hygiene. The booklet also answers the most frequently asked questions, thus helping the patient in adapting to using continence products.

### Accessories

Along with a selection of urine bags, Flexicare home packs also contain the following accessories:

1 x Pair of anti-slip latex-free lycra velcro leg straps

1 x FlexiHang 2 litre bag hanger

1 x User booklet

1 x Prescription card, with details of all products available and prescription order numbers

“Providing ongoing care in the community”



00-0040	1 x F4, 1 x F2 and 1 x 500ml short tube leg bag - continence care home pack	1
00-0042	4 x F2 and 1 x 500ml short tube leg bag - continence care home pack	1
00-0045	2 x F4L and 1 x 500ml short tube leg bag - continence care home pack	1
00-0046	2 x F4L and 1 x 500ml long tube leg bag - continence care home pack	1
00-0047	2 x F2S and 1 x 500ml short tube leg bag - continence care home pack	1
00-0048	2 x F2S and 1 x 500ml long tube leg bag - continence care home pack	1

### Headquarters

Flexicare Medical Limited  
Cynon Valley Business Park,  
Mountain Ash, Mid Glamorgan.  
CF45 4ER, United Kingdom  
☎ +44 (0)1443 474647  
☎ +44 (0)1443 474222  
✉ enquiries@flexicare.com

### Flexicare Middle East

Flexicare Middle East Amman  
Um Uthina Shatt Al Arab Street  
Nakho Business Center  
2nd Floor - Office No.5  
☎ +(962) 6 554 89 44  
☎ +(962) 6 554 89 45  
✉ musa.wer@flexicare.com

### Flexicare India

Flexicare Medical India Pvt.Ltd  
North Phase,  
Developed Plot No. 10A,  
T.V.K. Industrial Estate,  
Ekkattuthangal,  
Chennai-600 097  
✉ enquiries@flexicare.com

### Flexicare China

Flexicare Medical China Limited  
No. B-15, Xicheng Industrial Zone  
I, Hengli Town, Dongguan City,  
Guangdong Province  
China 523460  
✉ enquiries@flexicareasia.com

### Flexicare USA

Flexicare Inc  
15 Hammond, Suite 305 Irvine,  
CA 92618 USA  
☎ 949 916 3333 (Toll Free)  
☎ 949 916 3303  
✉ inquiries@flexicare.com

**Flexicare Online** visit our website for all the latest news, events, educational resources and information on our full product range.

[www.flexicare.com](http://www.flexicare.com)

## We strive to exceed our customer's expectations in Quality, Service and Value.

In addition to supplying high quality products, we have a strong commitment to provide a first class service and achieve the highest standards of professionalism and efficiency in all areas of our operations.

We encourage new and innovative ideas from our customers and Flexicare's experienced technical and manufacturing team will be pleased to discuss developmental projects as well as custom products for your specific needs.

Flexicare's Customer Care Department is able to respond promptly to your enquiries and provide up to date information regarding our product range or to simply discuss your orders and their status.

Our team of Territory Sales Executives, with nationwide coverage supported by our Customer Care Department at Flexicare's headquarters, will be pleased to arrange to visit you to discuss the comprehensive range of Flexicare products or any other requirements you may have.



Flexicare Medical Limited. Cynon Valley Business Park, Mountain Ash, Mid Glamorgan, CF45 4ER, UK.

 +44 (0)1443 474647

 +44 (0)1443 474222

 [enquiries@flexicare.com](mailto:enquiries@flexicare.com)

 [www.flexicare.com](http://www.flexicare.com)